

Lumena™ (Straight-Tip) Ureteral Access Sheath

Provides reliable access to the ureter while protecting surrounding tissues and enhancing procedural control.

Ultra-Smooth Navigation


The specialized hydrophilic coating significantly reduces surface friction, enabling easier insertion while minimizing tissue damage.

Outstanding Kink Resistance

Incorporates a stainless-steel reinforced layer and anti-kinking design, enhancing durability and stability.




Lumena™ (Straight-Tip) Ureteral Access Sheath

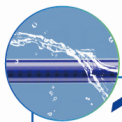


1 Ultra-Smooth Hydrophilic Coating

Enables effortless sheath placement through tight anatomical structures.

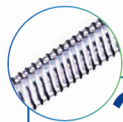


2 Sheath Tube with Reinforcing Coil

Enhances lumen integrity and prevents kinking, offering an excellent scope-to-sheath ratio.

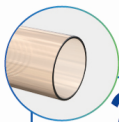

3 Ultra-Smooth PTFE Inner Wall

Ensures easier device insertion and removal, minimizing resistance.

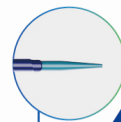

4 Soft Tapered Dilator Tip Design

Protects ureteral mucosa and surrounding tissues during insertion.

Product Specification

Model Specification	Effective Length L (mm)	Outer Diameter OD (mm)	Inner Diameter ID (mm)
UAS951135	350	3.85	3.20
UAS951140	400		
UAS951145	450		
UAS111335	350	4.45	3.85
UAS111340	400		
UAS111345	450		
UAS121435	350	4.70	4.05
UAS121440	400		
UAS121445	450		
UAS131535	350	5.00	4.40
UAS131540	400		
UAS131545	450		

*Full model specifications, options and reordering details can be found on our website or by reaching out to our support team.

Customizable Items

- ☐ Inner/outer packaging and accessories
- ☐ Product specifications (outer diameter, length, etc.)
- ☐ Bending length and angle
- ☐ Product appearance and batch number
- ☐ Catheter structure
- ☐ Others



CatheGenix Pte. Ltd.

+65 64655909
 info@CatheGenix.com
 www.cathegenix.com

2 Tukang Innovation Grove,
 #08-09/10 JTC Medtech Hub, Singapore 618305

