

NaviGlide™ Zebra GuideWire

Designed for precision, safety, and seamless navigation in urological procedures. It is utilized alongside an endoscope to direct catheters and insert instruments during urological surgery, with its superior visibility offering essential guidance for surgical procedures.

Enhanced Surface Performance

The PTFE exterior layer with hydrophilic coating achieves ultra-low friction, reducing the friction coefficient by 95%. This allows for significantly smoother navigation and improved device handling.


Tissue Protection Design

The 6 cm tip is coated with PU and a hydrophilic layer, providing effective protection for ureteral tissue and minimizing the risk of perforation.

NaviGlide™
SureTip

NaviGlide™
SoftTip


Before and after applying
hydrophilic coating on PU and PTFE.

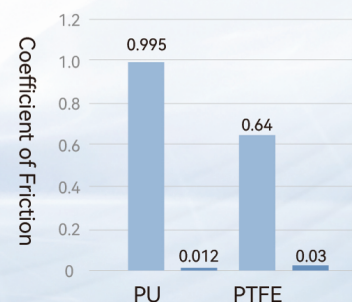

NaviGlide™ Zebra Guidewire

- 1 Anti-Twist Nitinol Core**
Provides kink resistance and 1:1 torque control for accurate device navigation.

- 2 Classic Blue-and-White Striped Design**
Enhances visual dynamics for clear positioning during procedures.

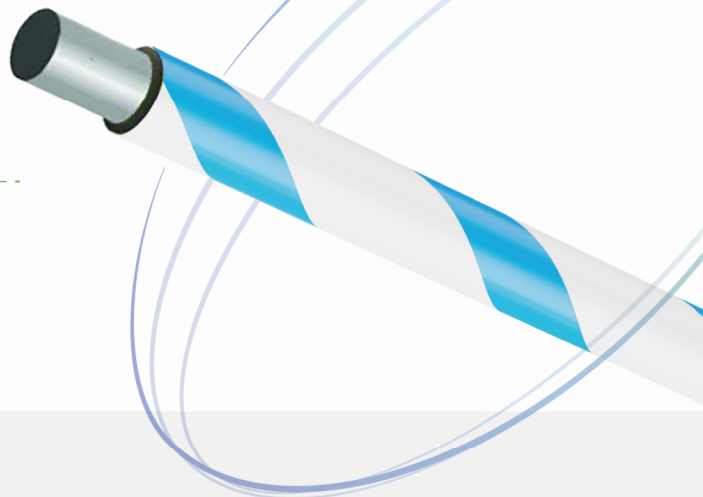

Product Specification

Guidewire Types	Diameter OD (inch)	Length L (cm)		
		150cm	260cm	450cm
SureTip™ Complete Striped Jacket	0.021	GWC021150	GWC021260	GWC021450
	0.025	GWC025150	GWC025260	GWC025450
	0.028	GWC028150	GWC028260	GWC028450
	0.032	GWC032150	GWC032260	GWC032450
	0.035	GWC035150	GWC035260	GWC035450
	0.038	GWC038150	GWC038260	GWC038450
SoftTip™ Striped Jacket with PU Tip	0.021	GWH021150	GWH021260	GWH021450
	0.025	GWH025150	GWH025260	GWH025450
	0.028	GWH028150	GWH028260	GWH028450
	0.032	GWH032150	GWH032260	GWH032450
	0.035	GWH035150	GWH035260	GWH035450
	0.038	GWH038150	GWH038260	GWH038450

*Full model specifications, options and reordering details can be found on our website or by reaching out to our support team.

Customizable Items

- ☐ Inner/outer packaging and accessories
- ☐ Guidewire specifications (outer diameter, length)
- ☐ Coating performance
- ☐ Physical properties (friction, pushability, kink resistance)
- ☐ Guidewire type
- ☐ Tip style and length
- ☐ Application scenario requirements, etc.



CatheGenix Pte. Ltd.

+65 64655909

info@CatheGenix.com

www.cathegenix.com

2 Tukang Innovation Grove,
#08-09/10 JTC Medtech Hub, Singapore 618305

